

Chengdu Besmile Biotechnology Co.,Ltd

Chengdu, China  
+86 28 85317108  
Info@cdbesmile.com

WWW.BSMDENTAL.COM

# Aconia<sup>®</sup>















## CATALOGUE








































WWW.BSMDENTAL.COM



# CAD/CAM MATERIAL

Aconia®	TT	SHT	ST	HT +	Standard & Master	Artist	WAX & PMMA & PEEK	Milling burs
Multilayer								
Preshade								
White								

# INDICATION

		Veneer	Inlay & Onlay	Reduced crown	Full contour crown		Coping	Full contour anterior bridge (3unit)	Full contour bridge (3unit)	Full contour bridge (≤7unit)	Full contour bridge (≤14unit)	Abutment
	<b>TT</b> (3D Multilayer & White)											
<b>HOT</b>	<b>SHT-ML</b> (3D Multilayer)											
	<b>SHT</b> (Preshade & White)											
	<b>ST</b> (Preshade & White)											
<b>UPGRADED</b>	<b>HT<sup>+</sup></b> (White)											


# TT-Multilayer

Multilayer

## TT-ML

Create the best smile with highest esthetics

- Most natural appearance
- Seamless transition of translucency, strength and shade
- Fast and easy processing
- Creatively maximized efficiency and esthetics

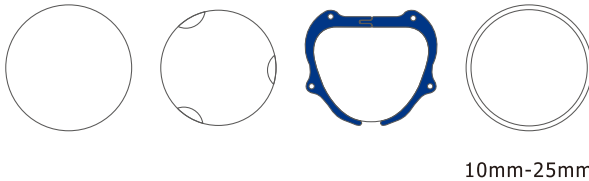

Indication					
Veneer	Inlay & Onlay	Reduced crown	Full contour crown	Coping	Full contour anterior bridge (3 unit)



- Seamless transition, smooth gradient
- Extremely high incisal translucency provides lifelike replication of tooth enamel
- Bionic tooth-growing effect created by Aconia Vitalization Technology

Technical data	
Flexural strength(3-point)	600-900 Mpa
Translucency	46-49%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>3/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>



# SHT-Multilayer

Multilayer

## SHT-ML

Create the best smile with highest flexibility

- All-in-one & One-for-all
- Seamless gradient in translucency, strength and shade
- Fast and easy processing
- Revolutionarily well-balanced combination of strength and translucency



Indication						
Inlay & Onlay	Reduced crown	Full contour crown	Coping	Full contour anterior bridge (3 unit)	Full contour posterior bridge (3 unit)	Full contour bridge (≤7 unit)



## ALL IN ONE

All in One

All technologies integrated to one **Aconia** masterpieces made in one Multilayer indications applied by one SHT-ML

One for All

One SHT-ML to fulfill All your needs  
Multilayer to replace all your inventories  
**Aconia** create all your smiles

Technical data	
Flexural strength(3-point)	900-1100 Mpa
Translucency	43-46%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>5/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>

# PRESHADE


## SHT-Preshade

Create the best smile with high efficiency

- Excellent strength combined with 46% translucency
- A wide range of indications for full contour restorations
- Reproduce Vita shades perfectly
- Efficient and ideal results



Technical data	
Flexural strength(3-point)	1000±100 Mpa
Translucency	46%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered Density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>5/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>

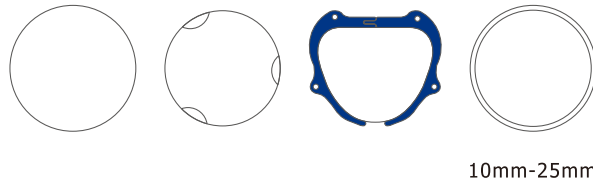

Indication						
Inlay & Onlay	Reduced crown	Full contour crown	Coping	Full contour anterior bridge (3 unit)	Full contour posterior bridge (3 unit)	Full contour bridge (≤7 unit)


## ST-Preshade

Create the best smile with speed

- Outstanding strength of 1250Mpa combined with ideal translucency
- A wide range of indications for restorations from coping to long-span bridge
- Reproduce Vita shades perfectly
- Stable and reproducible results



Technical data	
Flexural strength(3-point)	1250±100 Mpa
Translucency	43%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered Density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>5.5/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>

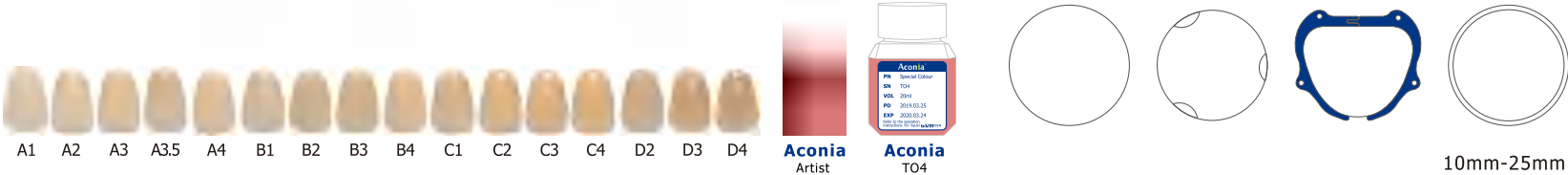

Indication						
Inlay & Onlay	Full contour crown	Coping	Full contour anterior bridge (3 unit)	Full contour posterior bridge (3 unit)	Full contour bridge (≤ 7unit)	Full contour bridge (≤14unit)

# WHITE



## TT

Create the bright smile: the professional solution

- Highest translucency up to 49%
- Artistic foundation for individualized restoration
- Esthetic alternative to lithium disilicate with double strength



Technical data	
Flexural strength(3-point)	> 700 Mpa
Translucency	49%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered Density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>3/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>

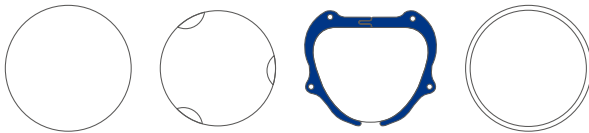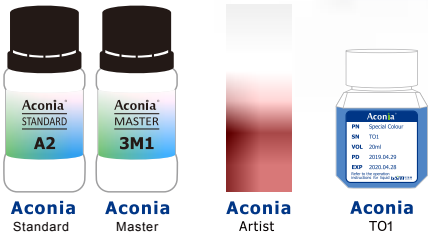

10mm-25mm

Indication					
Veneer	Inlay & Onlay	Reduced crown	Full contour crown	Coping	Full contour anterior bridge (3 unit)

# White



## SHT

Create the vivid smile: the cost effective solution

- Outstanding translucency of 46%
- Remarkable strength with reliable durability
- Wide indications for full contour
- Easy and fast coloring



Technical data	
Flexural strength(3-point)	1000±100 Mpa
Translucency	46%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered Density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>5/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>

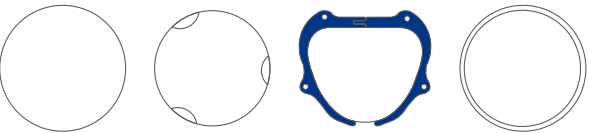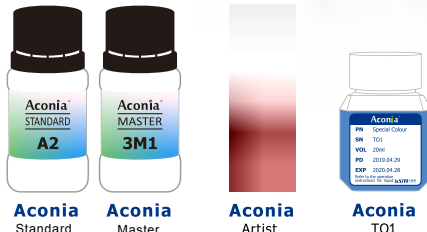

10mm-25mm

Indication						
Inlay & Onlay	Reduced crown	Full contour crown	Coping	Full contour anterior bridge (3 unit)	Full contour posterior bridge (3 unit)	Full contour bridge (≤7 unit)

# WHITE


## ST

Create the universal smile: the classic solution

- Strong with good millability
  - Attractive translucency with outstanding strength 1250Mpa
- Wide indications from coping to long-span bridge
  - Easy and fast coloring



Technical data	
Flexural strength(3-point)	1250±100 Mpa
Translucency	43%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered Density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>5.5/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>

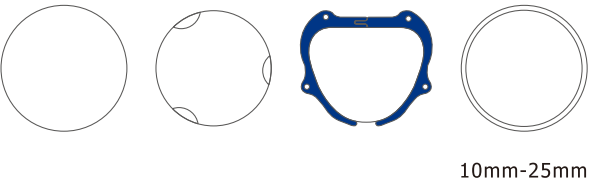

Indication						
Inlay & Onlay	Full contour crown	Coping	Full contour anterior bridge (3 unit)	Full contour posterior bridge (3 unit)	Full contour bridge (≤7 unit)	Full contour bridge (≤14 unit)

# White


## HT<sup>+</sup>

Create the unwavering smile: the economical solution

- The extraordinary strength of 1400Mpa grants a high level of process safety
- Easy veneering & individualizing
- Best option for coping & abutment



Technical data	
Flexural strength(3-point)	>1400 Mpa
Translucency	42%
Vickers-hardness HV10	1300±50
Density	>3 (g/cm³)
Sintered Density	>6.02 (g/cm³)
Chemical solubility	<50(µg/cm³)
Radioactivity	<0.1/Bq.g <sup>-1</sup>
Fracture toughness	>9/(Mpa.m <sup>1/2</sup> )
CTE	(10.5±0.5)*10 <sup>-6</sup> K <sup>-1</sup>

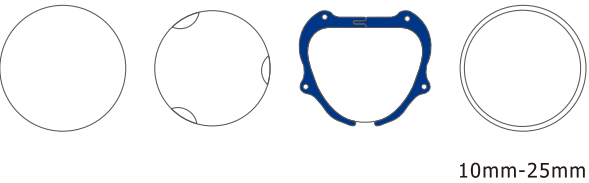

Indication							
Inlay & Onlay	Full contour crown	Coping	Full contour anterior bridge (3 unit)	Full contour posterior bridge (3 unit)	Full contour bridge (≤7 unit)	Full contour bridge (≤14 unit)	Abutment

# Standard & Master



50ml


## Standard

- Precisely matched to VITA\* 16 color system
- Well-suited for both dipping and brushing (paint-on) methods
- Fast-coloring & no color difference between pontic and neighboring crowns
- Ideal results applied on Aconia white blanks by Aconia Coloring Technology

The terms marked with \* are registered trademarks and/or brand names of the respective companies.



A1 A2 A3 A3.5 A4 B1 B2 B3 B4 C1 C2 C3 C4 D2 D3 D4

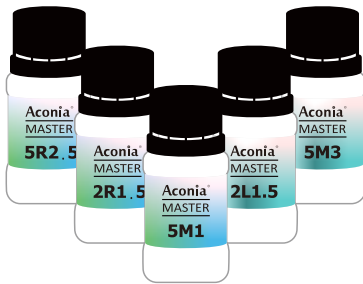

50ml

## Master

- Precisely matched to VITA\* 26 color system
- Well-suited for both dipping and brushing (paint-on) methods
- Fast-coloring & no color difference between pontic and neighboring crowns
- Ideal results applied on Aconia white blanks by Aconia Coloring Technology

The terms marked with \* are registered trademarks and/or brand names of the respective companies.



1M1 1M2 2L1.5 2L2.5 2M1 2M2 2M3 2R1.5 2R2.5 3L1.5 3L2.5 3M1 3M2 3M3 3R1.5 3R2.5 4L1.5 4L2.5 4M1 4M2 4M3 4R1.5 4R2.5 5M1 5M2 5M3


# Artist



## Artist

20ml

- Restore realistic, naturally lifelike appearance
- Reproduce rare and special colors
- Create esthetic art effect
- Enable individualized customization
- Start coloring process without preparing and mixing
- Deliver ideal results when applied on Aconia blanks

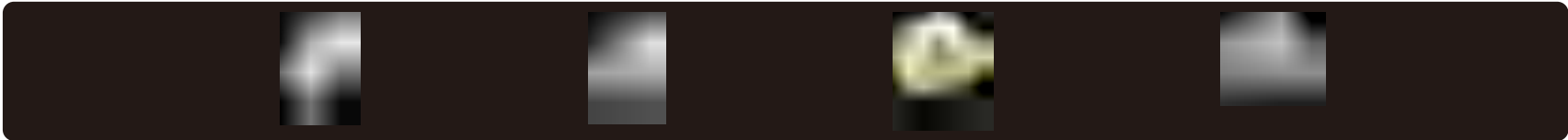

## Special colors

### Magic value



V1: Opaque



V2: Grey



V3: Blue

### Natural gingival



G1: Pink



G2: Red



G3: Violet

### Incisal translucency



TO2: Incisal Blue



TO3: Incisal Grey

### Art fissure



P1: Sunrise



P2: Sunset



P3: Olive



O1: Orange



O2: Brown




O3: Vanilla

# Wax



- Easy to mill
- High melting point
- Burn out completely without residue
- Compatible for making all male molds of oral tissue

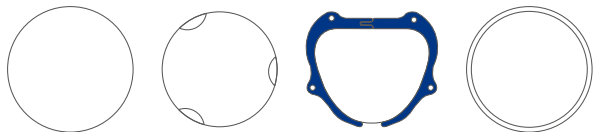


10mm-25mm

Technical data	
Material	Polymer
Color	Blue/White
Density	0.89-0.93g/m³
Drop melting point	106°C
Shore hardness	50-60 ShoreD

# PMMA



- Excellent abrasion resistance
- Excellent finishing bright surface & polishing performance
- Excellent long-term shade stability and esthetics
- To fabricate fully and partially long-term temporary crown and bridge

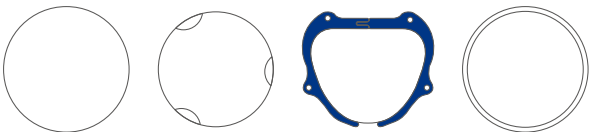

10mm-25mm

Technical data	
Material	PMMA 100%
Color	A0 A1 A2 A3 Pink Transparent
Density	1.19g/m³
Flexural strength	>125Mpa
Ash	<0.29%
Rate of contraction	<0.5%






# Milling burs






## CAD/CAM BURS: *Uncoated & Diamond*

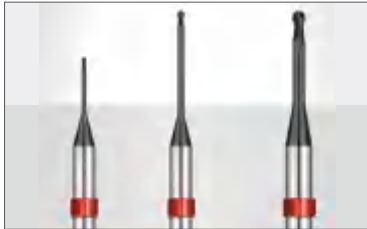- Excellent durability for the milling cutters
- Economical working- cost - benefit ratio
- With diamond coating for precise milling results and smooth surfaces.




## BSM BM-430D / BM-520D

	Head diameter	Diameter	Length
	2(R1.0)	4	22
	1(R0.5)	4	22
	0.6(R0.3)	4	12




## Roland

	Head diameter	Diameter	Length
	2(R1.0)	4	17
	1(R0.5)	4	15
	0.6(R0.3)	4	12




## Imes-icore 350i

	Head diameter	Diameter	Length
	2(R1.0)	6	17
	1(R0.5)	6	15
	0.6(R0.3)	6	12

## VHF K4/K5

	Head diameter	Diameter	Length
	2(R1.0)	3	17
	1(R0.5)	3	15
	0.6(R0.3)	3	12

## Imes-icore 240i/250i

	Head diameter	Diameter	Length
	2.5(R1.25)	3	17
	1(R0.5)	3	15
	0.6(R0.3)	3	12

# BM-430D



## 4 Axis milling machine

Besmile has ungraded 4 Axis BM-430D dental milling machine with a high-performance open system, which can meet the diverse needs of customers and ensure the long-term stability and accuracy.

- High speed data processing and analyzing NC system
- Compatible processing system
- Max.40,000rpm of the precise motorized spindle
- High-resolution digital servo-control system
- Automatic changing and calibrating of the 4 milling burs
- Visualized processing

## Technical data

L*W*H	600mm*470mm*650mm	Weight	70kg
Built-in burs	3 (2mm; 1mm; 0.6mm)	Linkage axis	4
Working Condition			
Voltage / power	220V//850W	Air pressure	>0.5MPa
Temp.	5°C~40°C	Motor type	Step-motor
Processing Category			
Millable material	Zirconium, PMMA, Wax, PEEK etc.	Millable Prosthesis	single crown, long bridge, inlay, onlay, veneering etc.
Milling accuracy	±0.01mm	Automatic tool-length measurement	Positive / √
Cooling system	Air-cooled spindle	Vacuuming system	Mute vacuuming*
Remote assistance	Positive / √	Dry milling	Positive / √
Touch-screen control	8" LED touch screen	Wet milling	Negative / X

\*match according to customer demand.

# BM-520D



## 5 Axis milling machine

The 5 Axis dental milling machine is developed independently by Besmile, which possesses high-precision mechanical structure and high-resolution control system, providing a brand new operation experience for users.

- Stable and reliable performance
- An extensive service-life
- Automatic 4 burs changing system
- Integrative positioning module
- Accuracy up to 5µm of servo-control system
- Real-time monitoring by the sound-sensing alarm
- Remote control & service system


## Technical data

L*W*H	680mm*570mm*780mm	Weight	150kg
Built-in burs	4 (2mm; 1mm; 0.6mm)	Linkage axis	5
Working Condition			
Voltage / power	220V/3.9Kw	Air pressure	>0.65MPa
Temp.	5°C~40°C	Motor type	Servo-motor
Processing Category			
Millable material	Zirconium, PMMA, Wax, PEEK etc.	Millable Prosthesis	Full denture, crowns, bridges, inlays, onlays, veneerings, etc.
Characteristic parameter			
Milling accuracy	±0.005mm	Automatic tool-length measurement	Positive / √
Cooling system	Full automatic circulating water-cooled spindle	Vacuuming system	Mute vacuuming*
Dry milling	Positive / √	Remote assistance	Positive / √
Wet milling	Negative / X	Touch-screen control	9.7" LED touch screen

\*match according to customer demand.



# CAD/CAM solution

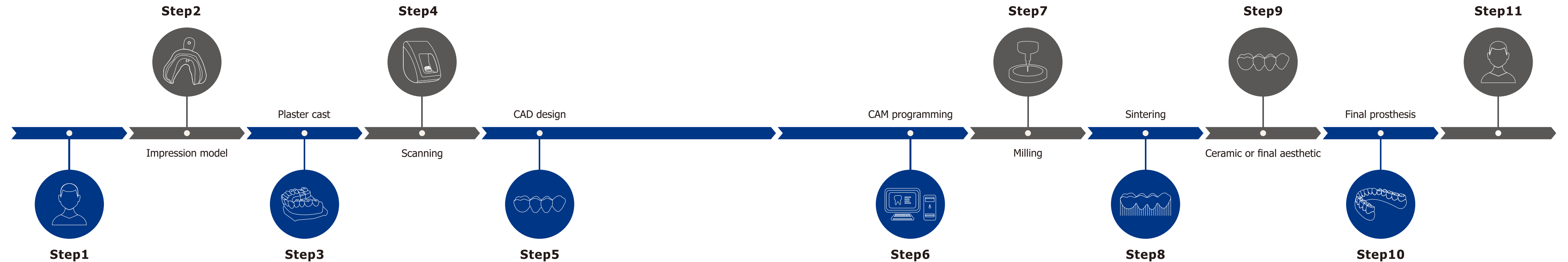

Besmile, dedicated in the dental zirconium blanks for dental **CAD/CAM solution**, is committed to delivering the beautiful smiles to dental patients worldwide by using its technology and reliable products together with efforts of dental technicians and dentists.

# Aconia®



Technology creates the best smile

